


roeko Surgitip ASPIRATOR TIPS

The product family for precise aspiration in surgery and endodontics



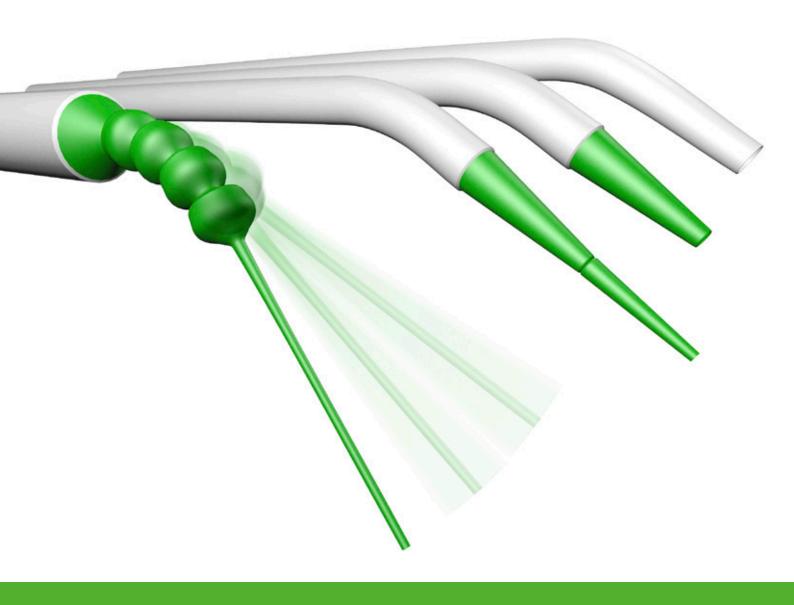

roeko

Surgitip Surgitip aspirator tips are proven aids in drying the oral field of treatment in minimally invasive surgery and in aspirating particles. The plastic aspirator tips allow extremely accurate aspiration as the shape and design have been specially developed for working under a microscope or loupe. The tips are ideal because of their hygienic handling and ease of use.


- The green tip provides excellent contrast for precise aspiration
- Open field of view due to the narrow shape and conical tip
- Hygienic and time-saving, being sterile, individually packaged and ready for immediate use
- Various aspirator diameters for a wide range of treatments



Surgitip-endo is an aspirating tip specially designed to dry root canals. The innovative design of the suction tip allows for a multi-directional flexibility, so it can easily be introduced into hard-to-reach root canals without having to bend the Canal Tip. Aspiration is highly effective at all angles and visibility is not constricted.

Surgitip-endo removes rinsing solutions and moisture quickly and effectively, shortening the drying time considerably. Less paper points are required. However, the use of paper points is recommended to check that the canal is thoroughly dry.

- Specially designed tip for multi-directional flexibility
- Effective aspiration and good canal accessibility
- Canal Tip ISO 60, inner diameter 0.35 mm
- Sterile, individually packaged and ready for use

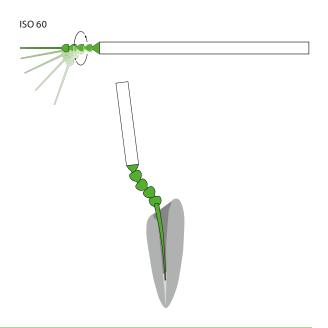

roeko

Surgitip-micro

With its 1.2 mm tip opening, the Surgitip-micro provides aspiration and precision during microsurgical procedures. The tip offers a good grip and is easy to manoeuvre, also making targeted aspiration possible in restricted dental operating environments. Surgitip-micro is sterile, individually packaged and ready for immediate use.



roeko Surgitip

Surgitip is used for surgical procedures that require greater aspiration capabilities. It is available with 2.5 mm and 4.8 mm tip openings. The 4.8 mm tip consists of the bent aspirator tube without the green tip. Surgitip is sterile, individually packaged and ready for immediate use.



roeкo

Aspirator tips


The non-sterile plastic aspirator tips are ideal for aspiration during non-surgical procedures. The length and shape correspond to those of the sterile 4.8 mm Surgitip. The aspirator tips are available in bulk packs with 100 pieces.



roeko

Double Adapter


The Double Adapter can be used to connect any ROEKO aspirator tip to standard aspirator systems with Ø 11 mm or Ø 16 mm tip openings. The Double Adapter is reusable and can be sterilised in the autoclave or chemiclave at up to 134 °C.

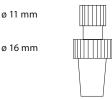

roeko

Concorde II

We recommend the Concorde II for universal aspiration. The plastic aspirator tip gently keeps the tongue and cheek away while aspirating large volumes. The tip made of high quality plastic is reusable and can be sterilised in the autoclave or chemiclave at up to 134 °C.

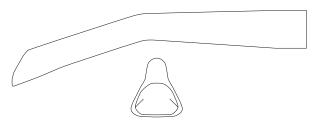

ORDER INFORMATION

	Tip opening ø mm	End opening ø mm	Contents	REF
Surgitip-endo	0.35 mm	6.5 mm	20 pcs*	6001 8280
Surgitip-micro	1.2 mm	6.5 mm	20 pcs * 100 pcs **	462 012 462 013
Surgitip	2.5 mm	6.5 mm	20 pcs ** 100 pcs **	462 025 462 026
	4.8 mm	6.5 mm	20 pcs * 100 pcs **	462 048 462 049
Aspirator tips	4.8 mm	6.5 mm	100 pcs *	461 048
Concorde		11 mm 16 mm	10 pcs *** 10 pcs	450 211 450 216
Adapter		11 mm 16 mm	1 1	463 025 463 011
Double Adapter		11 + 16 mm	1 Set	463 000

^{*}incl. 1 Double Adapter | **incl. 5 Double Adapters | ***incl. 1 Adapter





Coltène/Whaledent Ltd.

The President Suite Kendal House, Victoria Way Burgess Hill, West Sussex RH15 9NF / UK Tel +44 1444 235486 Fax +44 1444 870640 info.uk@coltene.com

Coltène/Whaledent

GmbH + Co. KG

Raiffeisenstraße 30 89129 Langenau / Germany Tel +49 7345 805 0 Fax +49 7345 805 201

info.de@coltene.com

Coltène/Whaledent Inc.

235 Ascot Parkway Cuyahoga Falls, Ohio 44223/USA Tel. + 1 330 916 8800 Fax + 1 330 916 7077 info.us@coltene.com

